

PHARMA SOLID KNOWLEDGE REPORT

Automated Process Development (APD)

Quality by Design (QbD) was introduced by the American FDA in the ICH Q8 (R2) "as a new paradigm with significant effects on current and future trends in pharmaceutical development" (2004). QbD describes good practice for pharmaceutical product development and introduces concepts of design space and flexible regulatory aspects. It is a systematic approach to development, emphasizes product and process understanding as well as process control.

Changes in formulation and manufacturing processes during development should be looked upon as opportunities to gain additional knowledge and further possibilities for improvements (from ICH Q8 (R2)).

WHAT IS QbD?

- Systematic approach to development
- Begins with predefined objectives
- Emphasizes product and process understanding and process control
- Based on sound science and quality risk management


BENEFITS OF QbD

- Enhanced process understanding
- Higher process capability
- Better product quality


Source: <https://www.fda.gov/media/85369/download>




CRITICAL QUALITY ATTRIBUTES AND QUALITY BY DESIGN

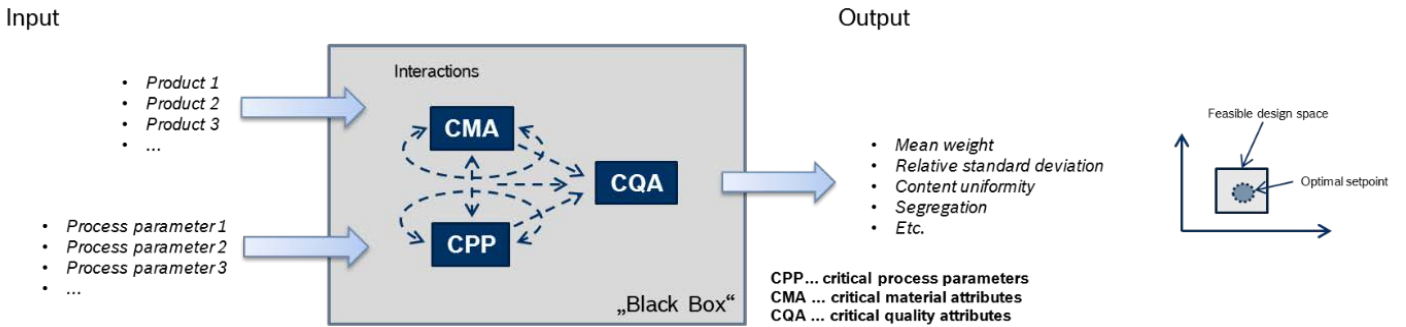


1. PLAN

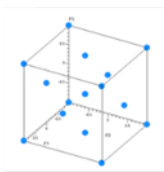


2. AUTOMATED EXECUTION

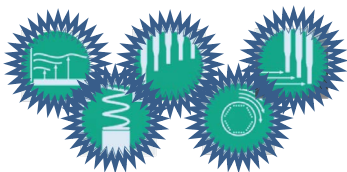

3. BUILD & ANALYZE PROCESS MODEL


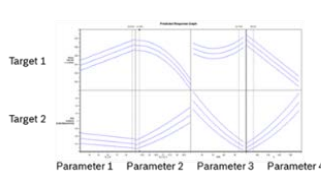


4. OPTIMIZE PROCESS & SOLVE PROBLEM

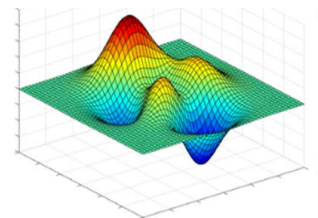


The combination of systematic planning (DoE) with automated test setup on a Syntegon GKF capsule filler is a unique and patented technology to perform many experiments in a short time. The automation allows you a broad screening to identify interactions between process parameters and to optimize yield, quality and robustness. The reproducibility of the tests and the quality of the data analysis is improved by an order of magnitude. Automation of all critical process parameters like power bed height, pneumatic spring rate or immersion depth combined with additional PAT-sensor technology and an automation control allow a fast, efficient and powder saving execution of systematically planned experiments.


Process Optimization

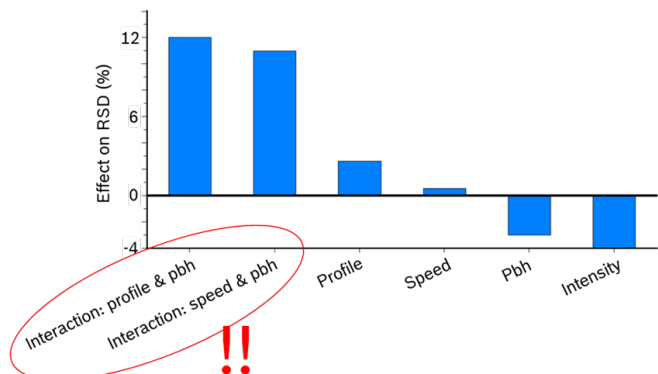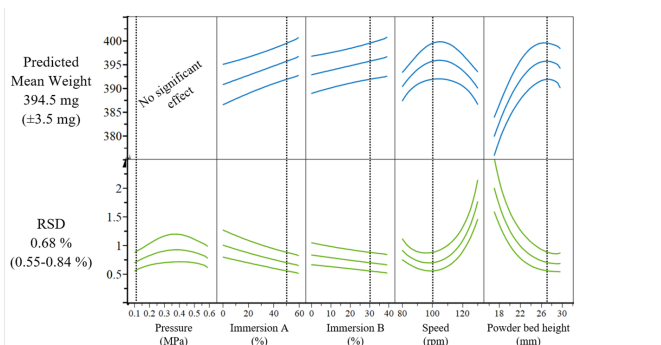

Formulation Development

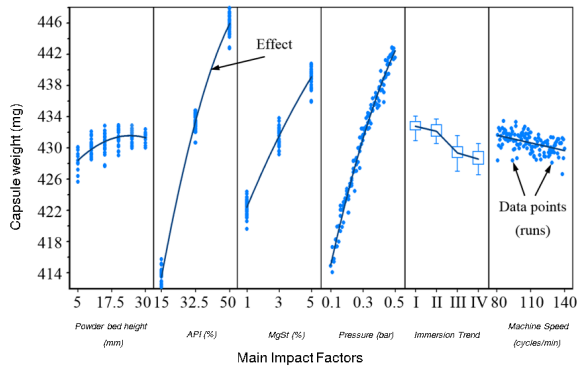

Identifying Interactions and CPP

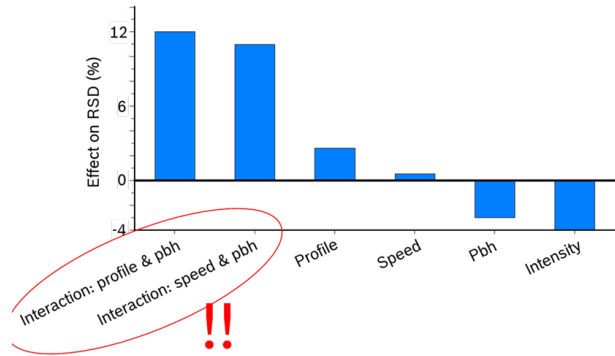


Excipient Investigation

Excipient	Type	Speed (rpm)	Pressure (MPa)	Powder bed (mm)	Height 1 (mm)	Height 2 (mm)	RSD (%)
Pearlit® Flash	Mannitol 1	120	0.08	3.1	5.0	3.1	1.1
Pearlit® 50C	Mannitol 2	120	0.18	14.8	3.2	3.0	1.2
Pearlit® 160C	Mannitol 3	121	0.21	15.0	3.2	5.0	1.3
Pearlit® 200SD	Mannitol 4	96	0.08	3.0	1.0	1.8	1.3
Manuogem EZ	Mannitol 5	129	0.08	14.5	1.0	2.8	1.5
Manuogem 2080	Mannitol 6	136	0.08	3.0	1.0	3.1	1.6
Paracel MD9	Mannitol 7	80	0.50	15.0	5.0	1.0	1.5
Paracel Xpress	Calcium lactate 1	120	0.35	15.0	5.0	3.2	0.7
Paracel DC	Calcium lactate 2	140	0.50	13.9	1.0	5.0	1.5
Paracel PP	Calcium lactate 3	120	0.17	15.0	5.0	2.3	1.4
Calcium lactate pentahydrate	Calcium lactate 4	127	0.35	15.0	3.1	1.0	1.3
Dicalcium phosphate	Calcium phosphate 1	140	0.08	4.2	1.0	3.3	1.7


Process parameter setting: low, middle, high

MORE THAN 100 APD TESTS PERFORMED!

- ❑ **RSD:** RDS of difficult sticky powder with parameter interaction improved from 9% to 2%
- ❑ **Weight:** 200mg transferred from capsule size 0 to capsule size 1 by improved tamping process parameters for a fluffy product
- ❑ Increased **yield** of good capsules from 94% to 99% for difficult products
- ❑ **Tool** for targeted **scale-up-process** and ability to handle batch variability by identifying CPP's
- ❑ Coming soon: APD for micro dosing of powders
- ❑ Coming soon: APD for TPR tablet compression

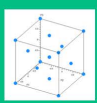



Design of Experiment (DoE)




Real machine experiments

Quality by Design with Automated Process Development (APD)

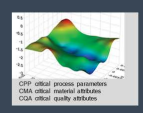

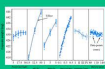

Integrated Measurement of the quality attributes



Automation & Automated testing

Conclusion

- ❑ The Automated Process Development (APD) is a unique Syntegon approach to real Quality by Design and optimum processes
- ❑ APD allows an efficient and fast determination of Critical Material Attributes, Critical Process Parameters and Critical Quality Attributes
- ❑ Enhanced process understanding and higher product capability leads to better product quality
- ❑ Simplified and controlled scale up process

Please contact us!

Please share your difficult products and processes with us, to find optimized solutions. Our "Engineering Pharmaceutical Service" team will be available with all our experience of over 50 years:

Dr. Thomas Brinz
Head of Department Engineering Pharmaceutical Service
Phone: +49(7151)14-2160
Thomas.Brinz@syntegon.com

Syntegon Technology GmbH

Stuttgarter Straße 130

71332 Waiblingen

Germany

Mail packaging-ph@syntegon.com

Web www.syntegon.com

facebook.com/Syntegon

twitter.com/Syntegon

youtube.com/Syntegon

linkedin.com/Syntegon